

1155 Out of Bed Monitor (radio)

1155 advantages verses traditional mats

There are inherent drawbacks to using floor mats to monitor out of bed residents, these are;

Mats

They are wired to a fixed point.
They tend to "creep" and cause broken wires.
They have a fixed service life of two years.
Residents can avoid stepping on them.
Two mats are needed to cover both sides of the bed.
Difficult to move for cleaning.
Specialist installers needed.

The benefits of using the 1155 are;

Radio, so no wires.
Portable, so easy to move to another room.
Monitors the whole area around the bed.
No cleaning issues.
Indefinite service life.
Low battery warnings
Out of range warnings
Tamper warning

A4721: 1155 Movement Monitor

The 1155 floor standing movement monitor is placed facing the bed of the person to be monitored. A single SET button is pressed on top of the PAM and the Carer leaves the room. Should the person get out of bed the PAM will transmit an "Out of Bed" signal to the 1133 panel and raise an alarm.

Features

- Easy to install and use.
- Rechargeable batteries last up to 4 weeks.
- Single button operation.
- Robust case
- Sealed to IP55
- Automatic low battery warning .
- 868Mhz radio frequencies allow longer range up to 80 metres indoor.
- Automatic out of range alarm.
- 30 second exit timer.



1142 Master Panel



1155 out of bed monitor



Key Points

- Simple plug and play installation
- Long battery life
- Bright OLED Display with integrated Sounder
- Sealed to IP55
- 868Mhz Radio allowing longer indoor range
- Automatic low battery warning
- Range of useful portable accessories


1155 Out of Bed Monitor (radio)

A4721: 1155 Movement Monitor

The 1155 floor standing movement monitor is placed facing the bed of the person to be monitored. A single SET button is pressed on top of the monitor and the Carer leaves the room. Should the person get out of bed the monitor will transmit an "Out of Bed" signal to the display panel and raise an alarm.

Features

- Easy to install and use.
- Rechargeable batteries last up to 4 weeks.
- Single button operation.
- Robust case
- Sealed to IP65
- Automatic low battery warning .
- 868Mhz radio frequencies allow longer range up to 80 metres indoor.
- Automatic out of range alarm.
- 30 second exit timer.

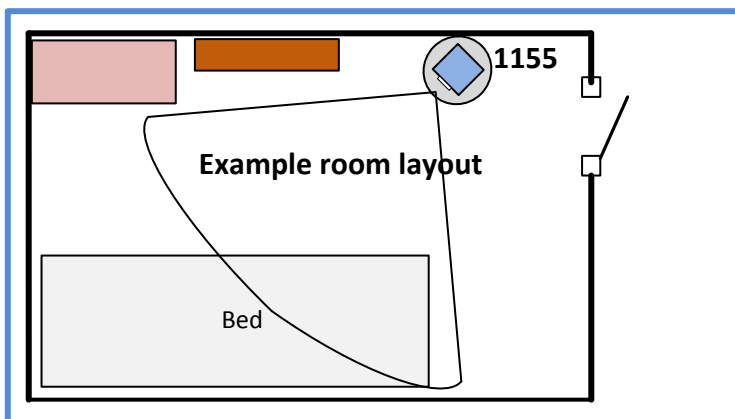

A6168: 1142 Radio Staff Receiver Display

The 1142 Master panel will accommodate up to 50 monitors. Operating by radio makes it easy to install and use. The 1142 when used with the 1155 floor standing movement monitor provides a simple to use system that is ideal for Carers & Nurses to use. The display will show how many 1155 monitors are in use and the signal strength as well as which monitors have a low battery.

For your peace of mind the 1142 master panel monitors the 1155 monitor radio signal and will set an alarm when any monitor is out of range.

- Receives alerts from up to 50 call units.
- A clear OLED display showing call location.
- High level integrated sounder.
- Multiple displays for larger premises.
- Power supply with backup battery.
- 868Mhz radio frequencies allow longer range up to 80 meters indoor.
- Portable master panels if required.

Key Points

- **Simple plug and play installation**
- **Long battery life**
- **Bright OLED Display with integrated Sounder**
- **Sealed to IP64**
- **868Mhz Radio allowing longer indoor range**
- **Automatic low battery warning**
- **Range of useful portable accessories**